
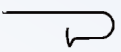
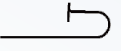
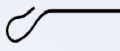
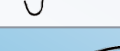
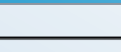


Product Specifications

Tip Configuration	Tip Shape	Size/Length(cm)	Tip Configuration	Tip Shape	Size/Length(cm)
	RH	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		HEADHUNTER 1	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	YAS	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		HEADHUNTER 3	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	RS	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		HEADHUNTER 5	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	RLG	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		HEADHUNTER 6	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	COBRA 1	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		NEWTON 1	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	COBRA 2	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		NEWTON 2	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	COBRA 3	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		NEWTON 3	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	Multifunctional tubeA 1	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		NEWTON 4	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	Multifunctional tubeB 1	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		NEWTON 5	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	SIMMONS 1	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		BENTSON 1	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	SIMMONS 2	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		BENTSON 2	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	SIMMONS 3	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		BENTSON 3	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	SIMMONS 4	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		MANI	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	Pulmonary artery	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		RENAL DOUBLE CURVE	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	STRAIGHT CATHETHERS	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		MIKAELSSON	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	RUC	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		SHEPHERD HOOK 0.8	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	REUTER	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		SHEPHERD HOOK 1.0	4F, 5F, 6F, 7F 60/65/80/90/100/110/125
	Single Curve	4F, 5F, 6F, 7F 60/65/80/90/100/110/125		VERTEBRAL	4F, 5F, 6F, 7F 60/65/80/90/100/110/125

Zhejiang Barty Medical Technology Co.,Ltd

Add:Room 101, 201, 301, Building 26, Yinhai Technology Innovation Center, Xiasha Street,
Qiantang District,310018 Hangzhou City, Zhejiang Province, China
Tel:0571-88373719
Fax:0571-88373716
Web: <https://www.barty-medical.com/>
Email: info@bartymedical.com

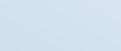

PENTY™
Angiography Catheter




PENTY™

Angiography Catheter




Wide Inner Lumen

Wide inner lumen, ensuring high efficiency flowing of contrast agent.

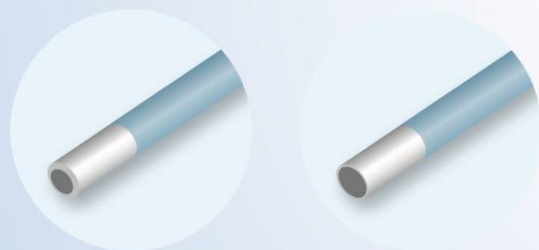

Stainless Steel Braided Net

Unique stainless steel braided net provides excellent 1:1 torsion control and kink resistance.

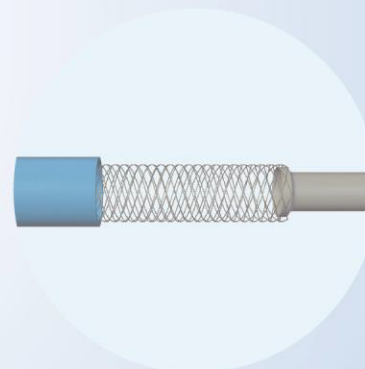

Smooth Outer Wall

Tubing outer wall polished to be smooth, allowing good propelling and trafficability.

Smooth Transition of Three-segment Design

Distal end: soft and mild, with excellent radiopacity;
Remote end: Optimize softness of the head, ensuring gentleness of operation and coaxality of blood vessel;
Body: Hard but submissive, ensuring an exact torsion transmission

Application

This product is used to inject radiographic media for peripheral vascular system during interventional operation.