


Product Specifications


Tip Configuration	Effective Length	5Fr	6Fr	7Fr	8Fr
 RDC	60cm	D-RDC-05060	D-RDC-06060	D-RDC-07060	D-RDC-08060
	65cm	D-RDC-05065	D-RDC-06065	D-RDC-07065	D-RDC-08065
	80cm	D-RDC-05080	D-RDC-06080	D-RDC-07080	D-RDC-08080
	90cm	D-RDC-05090	D-RDC-06090	D-RDC-07090	D-RDC-08090
	100cm	D-RDC-05100	D-RDC-06100	D-RDC-07100	D-RDC-08100
	110cm	D-RDC-05110	D-RDC-06110	D-RDC-07110	D-RDC-08110
	125cm	D-RDC-05125	D-RDC-06125	D-RDC-07125	D-RDC-08125
 MP	60cm	D-MP-05060	D-MP-06060	D-MP-07060	D-MP-08060
	65cm	D-MP-05065	D-MP-06065	D-MP-07065	D-MP-08065
	80cm	D-MP-05080	D-MP-06080	D-MP-07080	D-MP-08080
	90cm	D-MP-05090	D-MP-06090	D-MP-07090	D-MP-08090
	100cm	D-MP-05100	D-MP-06100	D-MP-07100	D-MP-08100
	110cm	D-MP-05110	D-MP-06110	D-MP-07110	D-MP-08110
	125cm	D-MP-05125	D-MP-06125	D-MP-07125	D-MP-08125

Zhejiang Barty Medical Technology Co.,Ltd

Add:Room 101, 201, 301, Building 26, Yinhai Technology Innovation Center, Xiasha Street,
Qiantang District,310018 Hangzhou City, Zhejiang Province, China
Tel:+86 0571-88373719
Fax:+86 0571-88373716
Web: <https://www.barty-medical.com/>
Email:info@bartymedical.com

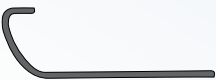

CONTY™
Guide Catheter

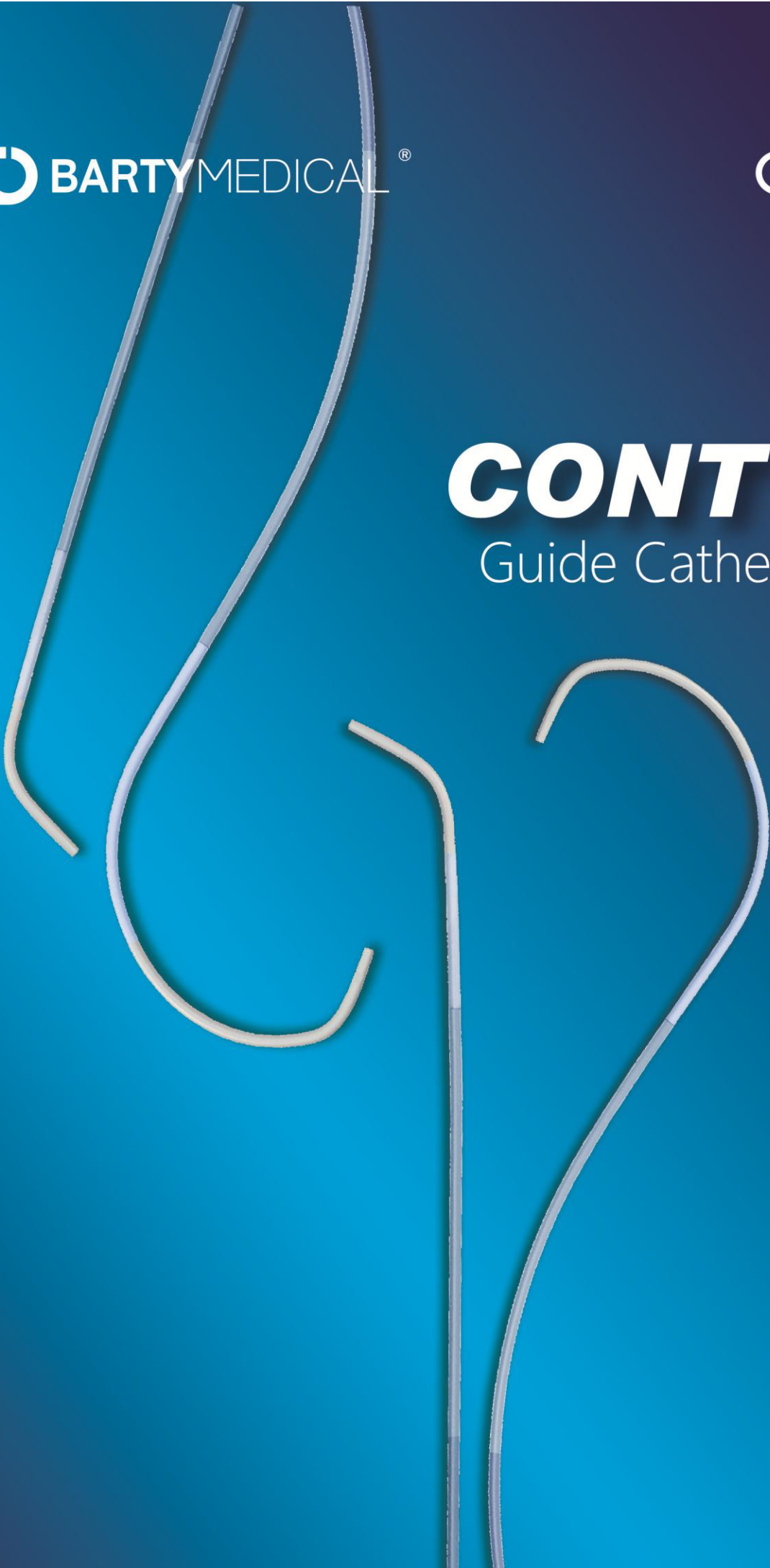




CONTY™

Guide Catheter


• Wide Inner Lumen

With wide inner lumen, it's easier for guidewire, balloon, stent, etc. to pass through.

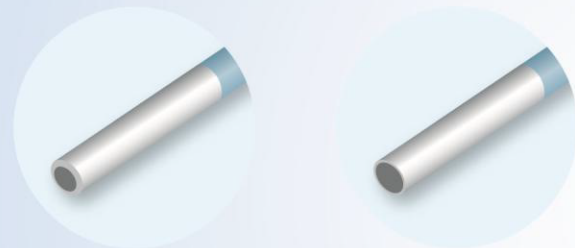

• PTFE Coating

PTFE coating over inner layer to reduce frictional resistance, providing a smooth delivery for device.

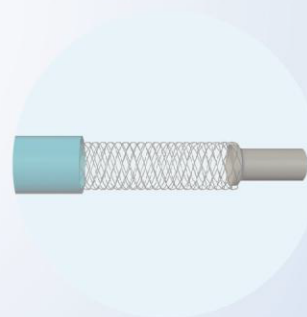

• Smooth Outer Wall

Tubing outer wall polished to be smooth, allowing good propelling and trafficability.

Application

This product provides a channel for balloon dilatation catheter, guidewire and other medical devices to insert during interventional operation, being used for peripheral vascular system.

Smooth Transition of Four-segment Design

Distal end: soft and mild, with excellent radiopacity;
Remote end: Optimize softness of the head, ensuring gentleness of operation and coaxality of blood vessel;
Proximal end: Nylon polymer material provides strong support
Body: Hard but submissive, ensuring an exact torsion transmission